

ADVANCED MANUAL INFLATE UA-705 Blood Pressure Monitor

Fast . Easy . Accurate .



A manual inflation digital blood pressure monitor with advanced features.


The Advanced Manual Inflate is an innovative monitor that incorporates digital blood pressure technology and advanced features at an affordable price. The monitor stores the last 30 readings in memory and also displays average readings. The average reading feature automatically calculates the “average” of the total readings stored in memory.

The Advanced Manual Inflate monitor also offers two innovative features, the Pressure Rating Indicator™ and Irregular Heartbeat feature. The Pressure Rating Indicator provides a snapshot of your blood pressure classification based on your measurements and eliminates the guesswork on what your blood pressure measurements mean. In addition, this monitor alerts the user of the presence of an Irregular Heartbeat and provides an accurate blood pressure and pulse rate measurements even when an Irregular Heartbeat occurs.

- Requires only one ‘AA’ battery (included)
- Available in medium or large cuff sizes (cuffs can also be purchased separately)
- Large-print instruction guide with illustrations and quick start card in English, French and Spanish
- Four week blood pressure logbook included

For questions about LifeSource products, call us toll free at 1-800-726-3364 and a specially trained representative will assist you, or visit us at www.LifeSourceOnline.com.

- One button operation
- 30 reading memory
- Pressure Rating Indicator™
- Irregular Heartbeat feature
- Displays average readings
- Large digital display
- Fast measurement
- Latex free

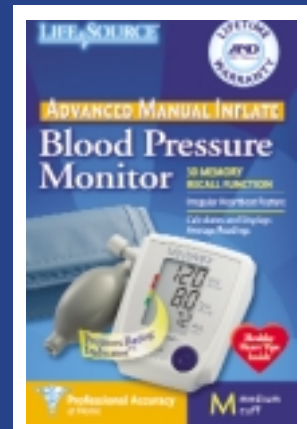

Advanced Manual Inflate UA-705 Blood Pressure Monitor

To obtain an accurate reading, it is important to use the correct cuff size. A cuff that is too small will result in a higher than expected reading and a cuff that is too large will result in a lower than expected reading. The Advanced Manual Inflate unit comes in two different cuff sizes to accommodate your customer's needs. This monitor is not compatible with the LifeSource small cuff.

	PART NUMBER	UPC NUMBER	CUFF SIZE
UNITED STATES	UA-705V	0-93764-60157-6	Medium 9.4 - 14.2" (24 - 36 cm)
	UA-705VL	0-93764-60158-3	Large 14.2 - 17.7" (36 - 45 cm)
CANADA	UA-705CN	0-93764-60159-0	Medium 9.4 - 14.2" (24 - 36 cm)
	UA-705CNL	0-93764-60160-6	Large 14.2 - 17.7" (36 - 45 cm)

Each unit comes with a trilingual instruction manual, quick start card, AA battery, four week blood pressure logbook and warranty card.

Accessories

PART #	DESCRIPTION
UA-700	Gray Latex Free Inflator Bulb
PA:4000735	Exhaust Valve
UA-280	Medium Blood Pressure Monitor Cuff: 9.4" - 14.2" (24-36cm)
UA-281	Large Blood Pressure Monitor Cuff: 14.2" - 17.7" (36-45cm)

Technical Specifications

Dimensions	3.2" (w) x 4.2" (d) x 2.2" (h) (81mm x 105mm x 54mm)
Weight	4.3 oz. (120g) without batteries
Display	Digital: Pressure - 15mm character height; Pulse - 10mm character height. Pressure and pulse displayed simultaneously.
Pressurization	Manual, by inflator bulb
Deflation	Hand-held push button
Type	Oscillometric
Accuracy	Pressure: ±3mmHg or 2%, whichever is greater/Pulse: ±5%
Pressure Range	20mmHg to 280mmHg/Pulse: 40 to 200 pulses/minute
Power Source	One type AA alkaline batteries (1.5 volt). Battery included.
Battery Life	Approx. 2,000 measurements
Operating Environ.	50°F to 104°F (10°C to 40°C)
Storage Environ.	14°F to 140°F (-10°C to 60°C)

Blood pressure measurements determined with this device are equivalent to those obtained by a trained observer using the cuff/stethoscope auscultation method within the limits prescribed by the American National Standards Institute (ANSI/AAMI SP10) for electronic sphygmomanometers. Manufactured by an ISO 9000 and ISO 13485 registered company. Engineered in Japan. Made in China.

Packaging and Shipping Specifications

Selling Unit (Box)	W x D x H (inches)	WT (lb.)	Individual Item Cube (cubic feet)
All UA-705 models	5 1/16 x 3 1/2 x 7 11/16	1.32	0.08
Unit Case (10 units)	W x D x H (inches)	WT (lb.)	Case Cube (cubic feet)
All UA-705 models	18 1/2 x 11 x 8 1/2	14.34	1.00

Ask your A&D Sales Representative about point-of-purchase floor and counter displays available for this blood pressure monitor unit.



1555 McCandless Drive
Milpitas, CA 95035
www.LifeSourceOnline.com

LIFE SOURCE®



Pressure Rating Indicator is a trademark of A&D Medical. Patent Pending. LifeSource is a registered trademark of A&D Medical.
©2005 A&D Medical. All rights reserved.



Ask your A&D Sales Representative about other LifeSource products.

- Blood Pressure Monitors
- Digital Thermometers
- Personal Scales
- Stethoscopes

For questions about LifeSource products, call us toll-free at 1-800-726-3364 or visit our web site at www.LifeSourceOnline.com

LITLS705 8/05