

QUICK RESPONSE

UA-787EJ

Blood Pressure Monitor

Fast. Easy. Accurate.



- 60 reading memory
- Easy-Fit™ cuff fits most arms
- Date/Time Stamp
- Three programmable alarms
- Pressure Rating Indicator™
- Irregular Heartbeat feature
- Displays average readings
- Large digital display
- Fast measurement

A monitor with high-end features that's easy-to-use.


The Quick Response monitor blends digital blood pressure technology and advanced features with simple functionality. This monitor offers memory, time and date feature and more with the simple push of a button. The UA-787EJ stores the last 60 readings in memory and also displays average readings. In addition, the Quick Response offers a time and date function, which indicates when a measurement was taken. Users can also program three alarms per day to remind them when to take a reading.

Included in the monitor is the Easy-Fit™ Cuff, an innovative cuff that instantly fits your arm. It has been specifically designed to fit most arms (9-17" (23-43 cm)) and provide the highest level of comfort and accurate readings.

The UA-787EJ also offers two innovative features, the Pressure Rating Indicator™ and Irregular Heartbeat feature. The Pressure Rating Indicator provides a snapshot of your blood pressure classification based on your measurements and eliminates the guesswork on what your blood pressure measurements mean. In addition, this monitor provides blood pressure and pulse rate measurements even when an Irregular Heartbeat occurs.

- Carrying pouch and AC Adapter included
- Large-print instruction guide with illustrations and quick start card in English, French and Spanish
- Four week blood pressure logbook included

For questions about LifeSource products, call us toll free at 1-800-726-3364 and a specially trained representative will assist you, or visit us at www.LifeSourceOnline.com.

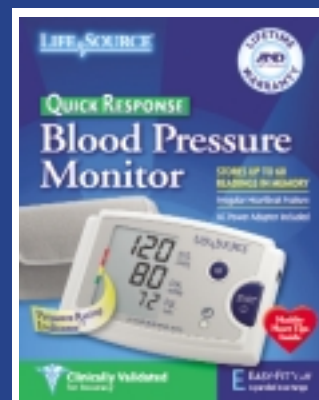

Quick Response UA-787EJ Blood Pressure Monitor

To obtain an accurate reading, it is important to use the correct cuff size. The Quick Response monitor includes the Easy-Fit™ Cuff. This cuff has been specifically designed to fit most arms (9-17" (23-43 cm)) and provide the highest level of comfort and accurate readings.

	PART NUMBER	UPC NUMBER	CUFF SIZE
UNITED STATES	UA-787EJ	0-93764-60154-5	9 - 17" (23 - 43 cm)
CANADA	UA-787EJCN	0-93764-60174-3	9 - 17" (23 - 43 cm)

Each unit comes with an AC Adapter, carrying pouch, trilingual instruction manual, quick start card, four week blood pressure logbook and warranty card.

Accessories

PART #	DESCRIPTION
TB:181A	AC Power Adapter

Technical Specifications

Dimensions	6.4" (w) x 4.4" (l) x 2.4" (h) (163mm x 112mm x 62mm)
Weight	12.3 oz. (350g) without batteries
Display	Digital: Pressure and pulse - 20mm character height Pressure and pulse displayed simultaneously.
Time/Date	12 hour format (AM/PM), year (2005-2054), month, and day with automatic adjustments for leap year
Pressurization	Automatic, using micropump
Deflation	Active exhaust valve (ECEVII)
Type	Oscillometric
Accuracy	Pressure: ±3mmHg or 2%, whichever is greater/Pulse: ±5%
Pressure Range	20mmHg to 280mmHg/Pulse: 40 to 200 pulses/minute
Power Source	Four type AA alkaline batteries (1.5 volt) (not included) or 120V AC Adapter.
Battery Life	Approx. 4 months with one daily measurement
Operating Environ.	50°F to 104°F (10°C to 40°C)
Storage Environ.	-4°F to 140°F (-20°C to 60°C)

This monitor is clinically validated according to the European Hypertension Society (EHS). Blood pressure measurements determined with this device are equivalent to those obtained by a trained observer using the cuff/stethoscope auscultation method within the limits prescribed by the American National Standards Institute (ANSI/AAMI SP10) for electronic sphygmomanometers. Manufactured by an ISO 9000 and ISO 13485 registered company. Engineered in Japan. Made in China.

Packaging and Shipping Specifications

Selling Unit (Box)	W x D x H (inches)	WT (lb.)	Individual Item Cube (cubic feet)
UA-787EJ	6 x 5 ³ / ₄ x 7 ¹ / ₂	2.38	0.15
Unit Case (10 units)	W x D x H (inches)	WT (lb.)	Case Cube (cubic feet)
UA-787EJ	32 x 11 ⁷ / ₈ x 11	28.40	2.62

Ask your A&D Sales Representative about point-of-purchase floor and counter displays available for this blood pressure monitor unit.



1555 McCandless Drive
Milpitas, CA 95035
www.LifeSourceOnline.com

LIFE SOURCE®



Ask your A&D Sales Representative about other LifeSource products.

- Blood Pressure Monitors
- Digital Thermometers
- Personal Scales
- Stethoscopes

For questions about LifeSource products, call us toll-free at 1-800-726-3364 or visit our web site at www.LifeSourceOnline.com