

AUTOMATIC WRIST UB-328 Blood Pressure Monitor

Fast. Easy. Accurate.




An automatic wrist monitor that is easy to use.

A&D Medical has applied its state-of-the-art measurement technology to monitoring blood pressure at the wrist. The Automatic Wrist monitor is activated with the push of a single button. The large display simultaneously shows systolic, diastolic and pulse readings.

This monitor stores the last 30 readings in memory and displays average readings. The average reading feature automatically calculates the “average” of all readings stored in memory. In addition, this monitor alerts the user of the presence of an Irregular Heartbeat and provides an accurate blood pressure and pulse rate measurements even when an Irregular Heartbeat occurs.

- Large-print instruction guide with illustrations and quick start card in English, French and Spanish
- Four week blood pressure logbook included
- Carrying case included

- Professional accuracy
- 30 reading memory
- Irregular Heartbeat feature
- Displays average readings
- Fast measurement
- Ideal for travel

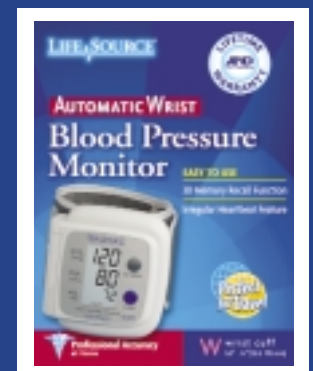

Automatic Wrist UB-328 Blood Pressure Monitor

The Automatic Wrist unit straps on quickly and comfortably to one's wrist to provide an accurate reading. This compact and portable unit gives fast measurements and provides convenience to those looking for an easy way to monitor their blood pressure.

	PART NUMBER	UPC NUMBER	CUFF SIZE
UNITED STATES	UB-328	0-93764-60036-4	5.3" - 7.7" (13.5-19.5cm)
CANADA	UB-328CN	0-93764-60137-1	5.3" - 7.7" (13.5-19.5cm)

Each unit comes with a trilingual instruction manual, quick start card, carrying case, four week blood pressure logbook and warranty card.

Technical Specifications

Dimensions	3.1 (l) x 2.9 (w) x 2.3" (h) (79 x 73 x 58 mm)
Dimensions with Case	4.0 (l) x 3.0 (w) x 3.3" (h) (102 x 76 x 83 mm)
Weight	4.2 oz. (120 g) without batteries
Display	Digital: 10/7 mm character height Pressure and pulse displayed simultaneously
Memory	30 readings
Pressurization	Automatic, using micropump
Deflation	Constant air release-valve system
Type	Oscillometric
Accuracy	Pressure: ± 3 mmHg or 2%, whichever is greater; Pulse: $\pm 5\%$
Pressure Range	20mmHg to 280 mmHg; Pulse: 40 to 200 pulses/minute
Power Source	Two type AAA alkaline batteries (1.5 volt). Batteries not included.
Battery Life	Approx. 400 measurements
Operating Environ.	50° to 104°F (10° to 40°C)
Storage Environ.	14° to 140°F (-10° to 60°C)

Blood pressure measurements determined with this device are equivalent to those obtained by a trained observer using the cuff/stethoscope auscultation method within the limits prescribed by the American National Standards Institute (ANSI/AAMI SP10) for electronic sphygmomanometers. Manufactured by an ISO 9001 and ISO 13485 registered company. Engineered in Japan. Made in China.

Packaging and Shipping Specifications

Individual Unit (Box)	W x D x H (inches)	WT (lb.)	Individual Item Cube (cubic feet)
All UB-328 models	4 ^{9/16} x 3 ^{7/16} x 6 ^{1/16}	0.65	0.06
Unit Case (10 units)	W x D x H (inches)	WT (lb.)	Case Cube (cubic feet)
All UB-328 models	18 x 10 x 7	7.75	0.73

Ask your A&D Sales Representative about point-of-purchase floor and counter displays available for this blood pressure monitor unit.



Ask your A&D Sales Representative about other LifeSource products.

- Blood Pressure Monitors
- Digital Thermometers
- Personal Scales
- Stethoscopes



1555 McCandless Drive
Milpitas, CA 95035
www.LifeSourceOnline.com

LIFE SOURCE®



LifeSource is a registered trademark of A&D Medical. ©2005 A&D Medical. All rights reserved.

LIFE:LS328 8/05

For questions about LifeSource products, call us toll-free at 1-800-726-3364 or visit our web site at www.LifeSourceOnline.com